WPB900

ELECTRIC BREAST PUMP

TIRE-LAIT ELECTRIQUE

ELEKTRİKLİ GÖĞÜS POMPASI

USER MANUAL **EN**

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KULLANIM KILAVUZU **TR**

Contents

Instruction manual	3
Important information	3
General information	3
Breast massage preparation.....	3
Names of parts	5
Dismantle directions	5
Assembly	6
Installation of exchanger locker.....	7
Power supply	8
Directions of breast pump using	8
Important	9
Cleaning the Pump	9
Troubleshooting	9
Storing breast milk.....	9
Warning	10
Repairs.....	10
Limited Warranty	10
Technical Specifications.....	10

Instruction manual

Dear breastfeeding mother,

In order to use this Electric Breast Pump effectively and safely, please read this manual carefully, and keep it in a safe place.

Important information

- *Please consult your doctor, lactation consultant, midwife or other breast feeding specialist to identify whether you are a suitable candidate for breastfeeding using breast pump.*
- *Never immerse the motor unit and AC adapter in water under any circumstances. intelligent*
- *Use only the original attachments supplied by the manufacturer, as only this will ensure proper functioning of the pump.*
- *If the pump does not function properly, any repairs must be carried out by professional staff*
- *Otherwise all guarantee claims are invalid.*
- *If the breast milk is for a premature baby, be sure to follow your doctor's instructions.*
- *Before using the breast pump please ensure that all components are clean and have been sterilised.*

General information

It is advisable to find a quiet, comfortable place where you can pump breast milk in a relaxed atmosphere without being disturbed.

Don't be surprised if the milk takes a little longer to start flowing properly at first. This is nothing unusual. In fact it is quite normal, since as a rule:

Breast milk production adjusts to suit daily requirements.

The daily amount required by infants can vary considerably. For this reason mothers are often recommended to drink as much as possible at first, in order to stimulate milk production.

You should pump milk first from one breast, then from the other. However, pumping should not exceed 7 minutes for each breast.

Breast massage preparation

Some Mothers may have difficulty in expressing milk either through using a manual breast pump or an Electric Pump. Using hot compresses with the massage method can assist in opening the mammary glands contributing to the flow of breast milk.

A hot compression can be achieved as follows:

- *Immerse a clean towel in warm water (60-70°C).*
- *Carefully squeeze out the excess water out of the towel.*

- Place the towel flat over your breasts until the towel cools down after a few minutes.
- Repeat steps 1-3 up to 4 times as required.

- After the hot compression has been completed, massaging the breasts may assist with opening the mammary glands.
- Start by placing your thumb and index finger around the base of the nipple and gently massage up and down the nipple about 5-6 times for each nipple.


- Move to the breast area and gently massage in a circular movement using four fingers flat and your thumb as shown in the picture. Stop after 5-6 rotations on each breast.

Note: If you experience pain during massaging or hot compression please stop immediately and seek medical advice.

Names of parts

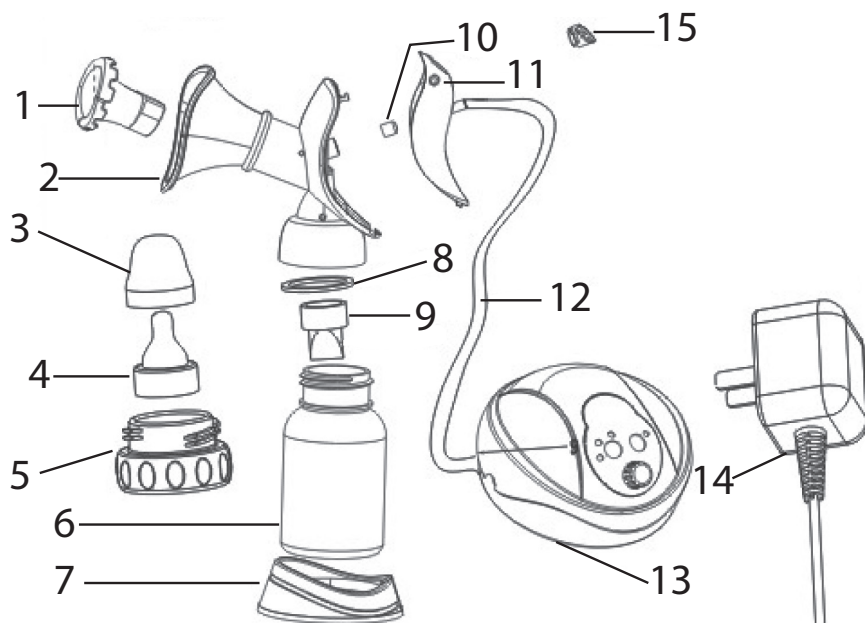- | | |
|------------------------|---------------------------|
| 1. Nipple protector | silicone gel/resist 110°C |
| 2. Breast shield cup | tritan/resist 110°C |
| 3. Teat cap | PP/resist 110°C |
| 4. Teat | silicone gel/resist 110°C |
| 5. Bottle adaptor | PP/resist 110°C |
| 6. Bottle | PP/resist 110°C |
| 7. Bottle stand | PP/resist 110°C |
| 8. Seal ring | silicone gel/resist 110°C |
| 9. Flow funnel | silicone gel/resist 110°C |
| 10. Seal plug | silicone gel/resist 110°C |
| 11. Electric connector | PP/resist 110°C |
| 12. Connect pipe | silicone gel/resist 110°C |
| 13. Motor unit | no sterilize |
| 14. AC adapter | no sterilize |
| 15. Locker | ABS/resist 110°C |

Dismantle directions

Before first use of this breast pump, please dismantle and sterilize all parts (except for motor unit and AC adapter):

- Turn this locker (15) anticlockwise to the maximum left position (figure 1).
- Hold breast shield cup (2) by left hand, dig the edge of middle part of this connector (11) by middle finger and index finger of right hand (figure 2).

- *Separate it from the breast pump shield (figure 3).*
- *Take out the nipple protector(1)*
- *Take out the seal plug(10)*
- *Dismantle bottle, teat and teat cap.*
- *Sterilize nipple protector, breast shield cup, teat cap, teat, bottle, bottle adapter, bottle stand.*

Figure 1


Figure 2


Figure 3

Assembly


- *Please ensure your hands are clean before touching the sterilised components. Insert the seal ring (8) and flow funnel (9) into the breast shield cup (2) tightly.*

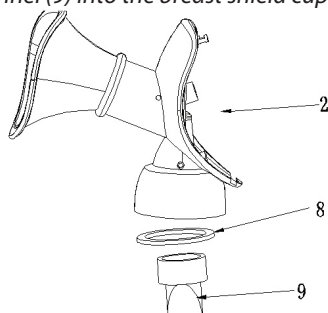
- *Insert the nipple protector (1) into the breast shield cup.*

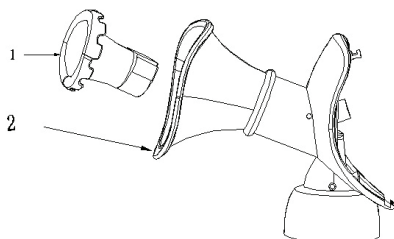- *Insert the flat end of seal plug (10) into the hole at back of breast shield cup (2), press it tightly, make sure there is no air gaps:*

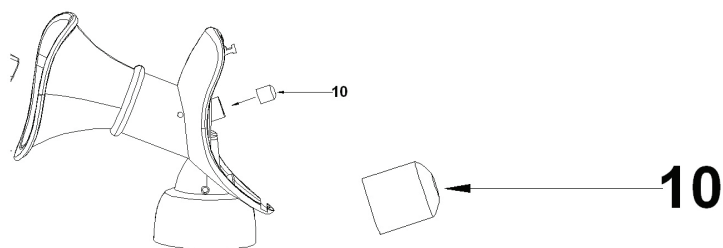- Turn locker (15) anticlockwise to the maximum left position (figure 4):
- Insert the locking pillar at back of the breast shield cup into the hole of exchange locker (figure 5). Next, push the middle of electric connector against the breast shield cup (figure 6), Next, press lower end of the electric connector and close those two items (figure 7).
- Turn the exchanger locker clockwise to the maximum right position (figure 8), lock two parts closely.
- Screw the breast shield cup onto the milk bottle and seal tightly.

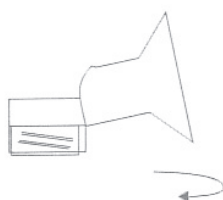- Connect this motor unit (13) and electric connector (11) with this silicone connect pipe (12)

Figure 4

Figure 5

Figure 6

Figure 7

Figure 8

Installation of exchanger locker

In case the locker (15) is squeezed out in accident, please install it as pictures follows:

- Hold the high side of the locker (figure 9)
- Aim the low side of locker at the gap in the round hole in the electric connector (figure 10).
- Press the exchange locker into the hole steadily. There will be crack sound (figure 11).
- Rotate the locker anticlockwise and clockwise to check whether or not it is installed well.

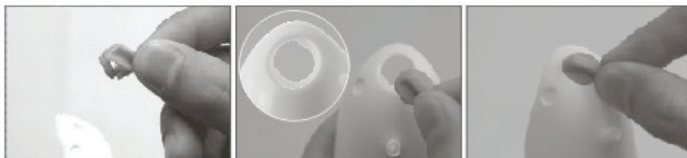
Figure 9

Figure 10

Figure 11

Power supply

- The electric breast pump can be run either with batteries or, using the adapter.
- For batteries operation you need four 1.5 V batteries, type AA. Make sure the batteries are the right way round.

Directions of breast pump using

- Rotate the suction adjuster dial to the maximum left position (figure 12).
 - Place the breast into the breast shield cup (figure 13), **with the nipple centered and sealed to the nipple protector (3)** (figure 14), to create a vacuum. Please ensure there are no air gaps between the nipple and the nipple protector, as well the breast and the breast shield cup.
 - Plug the unit into power or insert batteries and turn the switch in the motor unit to 'on' to start the motor.
 - Turn the dial at the front of the motor unit between 'L' and 'H' to get the proper suction (figure 12).
 - Press the 'Mode' button to get the most proper rhythm (figure 15).
- ◇ The pre-set mode is "Gentle". The suck lasts for 2 seconds. The mode light is blue colour.
 - ◇ Press one time, the mode is "Normal". The suck lasts for 3 seconds. The mode light is confused of blue and red.
 - ◇ Press two times, the mode is "Express", the suck lasts for 5 seconds. The mode light is red
 - ◇ Vaccum release lasts for 3 seconds for every mode.

Figure 12

Figure 13

Figure 14

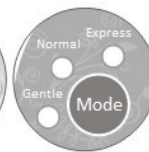
Figure 15

Important

- ***Make sure the nipple is centered in the nipple protector and sealed. Otherwise, there is risk that the milk enter into the Motor, which might damage the motor unit.***
- *During suck cycle, if you feel the suction is too strong, press “Mode” button to release vacuum. Move away the breast shield cup, adjust the suction to be proper.*
- *Always hold the breast pump in an upright position when using to ensure liquids do not enter into the motor unit.*

Cleaning the Pump

- *Dismantle all components and thoroughly wash (except for the motor unit and AC adaptor) using clean warm soapy water. The unit is also dishwasher compatible. Please ensure all components are placed in the top tray.*
- *To clean the motor unit simple wipe with a moist soapy cloth. Do not immerse the unit in water.*

Troubleshooting

Motor unit does not work

- *If using batteries please check the batteries are inserted correctly and are the right way around.*
- *Try replacing the batteries with super heavy duty alkaline or lithium batteries.*
- *Attempt to use with the adaptor.*

There is no suction

- *Turn the suction dial at the back of the unit to ‘High’.*
- *Check the battery strength (the batteries may need replacing).*
- *Ensure all components are fastened tightly and there are no air gaps.*
- *Ensure the breast is positioned in the breast shield cup so there are no air gaps.*
- *Ensure the flow funnel is placed and fixed firmly, no breakage on the soft edge*

I am not expressing milk

Some mothers may experience difficulty in expressing milk. Please first try the hot compression and massage techniques advised earlier in the manual. Also visit your breast feeding specialist to determine if you can breastfeed using an Electric Pump.

Storing breast milk

If the breast milk is not used immediately, firmly close the full bottle and keep it in the refrigerator.

If the breast milk is to be kept for more than 24 hours, it should ideally be deep-frozen (not in the pc bottle). We suggest you to use Weewell Breast Milk Bags.

Warning

Closed bottles filled with breast milk must not be placed in a microwave oven or deep-freeze. The strong heating/cooling could cause the bottle to explode. Uneven heating can result in scalding.

Repairs

This electrical appliance complies with relevant safety regulations. Repairs to electrical appliances must only be carried out by specialist personnel. Incorrect repairs can cause considerable danger to the user.

Limited Warranty

For a period of one year from the date of original purchase, Weewell warrants this product (excluding kit) against defects in material or workmanship. At Weewell's option, Weewell may repair or replace that is determined to be defective by Weewell and determined to be covered by this warranty. For warranty service, please return the product, along with dated sales receipt or other proof of the purchase date, to Weewell or to a local distributor. This warranty does not cover damage caused by accident, misuse, abuse, improper maintenance, unauthorized modification, or connection to an improper power supply. A charge will be made for repair of such damage. This warranty excludes all incidental or consequential damages and any liability other than what is stated above. Weewell makes no other warranties regarding this product (including warranties as to merchantability and fitness for a particular purpose) either express or implied.

Technical Specifications

Adaptor	Input: AC110~230V, 50~60Hz
	Output: DC5V 500~600mA
Batteries	4 x 1.5V AA Alkaline dry battery
Weight:	0.8kg
Operation Temperature:	>5°C

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Weewell Int'l Co., Ltd.

ADDRESS: 8F., No.45, Alley. 6, Lane. 182, Sec. 4, Chenggong Rd., Neihu Dist.,
Taipei City 114, Taiwan (R.O.C.)

Declaration of Conformity **Directive 99/5/EC (R&TTED)**

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We declare on our sole responsibility, that the following products:

Kind of equipments: Electric Breast Pump

Type-designation: Weewell-WPB900

are in compliance with the essential requirements of §3 of the R&TTED.

-Health and safety requirements pursuant to §3(1)a: (73/23/EEC)
Applied Standard(s) or other means of providing conformity:

EN 55014-1: 2006 +A1: 2009 +A2: 2011

EN 55014-2: 1997 +A1: 2001 +A2: 2008

EN 61000-3-2: 2006 +A1 +A2: 2009

EN 61000-3-3: 2008

-Protecting requirements concerning EMC §3(1)b: (2004/108/EEC)

We also hereby to confirm that the products, **WPB900** "Electric Breast Pump" is fully
complied with RoHS requirements.

Mr. Murphy Chen

Weewell Int'l Co., Ltd.

Weewell Int'l Co., Ltd.

ADRES: 8F., No.45, Ally. 6, Lane. 182, Sec. 4, Chenggong Rd., Neihu Dist.,
Taipei City 114, Taiwan (R.O.C.)

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Guarantee and After sales service

Guarantee conditions

This Weewell product is guaranteed, parts and labor, against all defects of conformity and latent defects for a period of two years from the initial purchase date of this product. This guarantee is applicable on presentation of a proof of purchase of the good.

This obligation of guarantee is committing Weewell or its legal representative only if the observed defect is not resulting from an irregular usage of the product or from a fall.

The repairs on products under guarantee can only be done by the Weewell after sales service or an authorized repairer. Other choice will imply the cancellation of the benefit of the legal two-year guarantee.

In case of problem

If the product is faulty, you will have to return it to your point of purchase, or you may contact Weewell after sales service by mail at: sav@weewell.com to find the nearest authorized repairer. Weewell technical assistance service in France can be contacted at: 33 (0) 800 101 291.

After sales service

Beyond the legal guarantee period, it is preferable to have any repair done through the Weewell after sales service or any authorized repairer, which addresses you will find in sending a mail to sav@weewell.com or in visiting the site www.weewell.com.

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Weewell Co., Ltd.